



PRESEPT™

EFFERVESCENT DISINFECTANT TABLETS



Hospital-grade disinfection validated to the latest **EU standards**


✓ Validated to 7 EN Standards

✓ Tested against 13 microorganisms

✓ EU BPR-compliant

✓ Suitable for critical healthcare environments

PRESEPT™ Effervescent Disinfectant Tablets are fully compliant with the EU Biocidal Products Regulation (EU BPR)¹ and validated to the latest EN standards. Their efficacy has been demonstrated across the defined EU resistance hierarchy, confirming performance against progressively more resilient pathogens relevant to healthcare settings. The following hierarchy reflects the standardized EU efficacy framework used to validate hospital disinfectants - supporting regulatory robustness and confidence in daily clinical use.

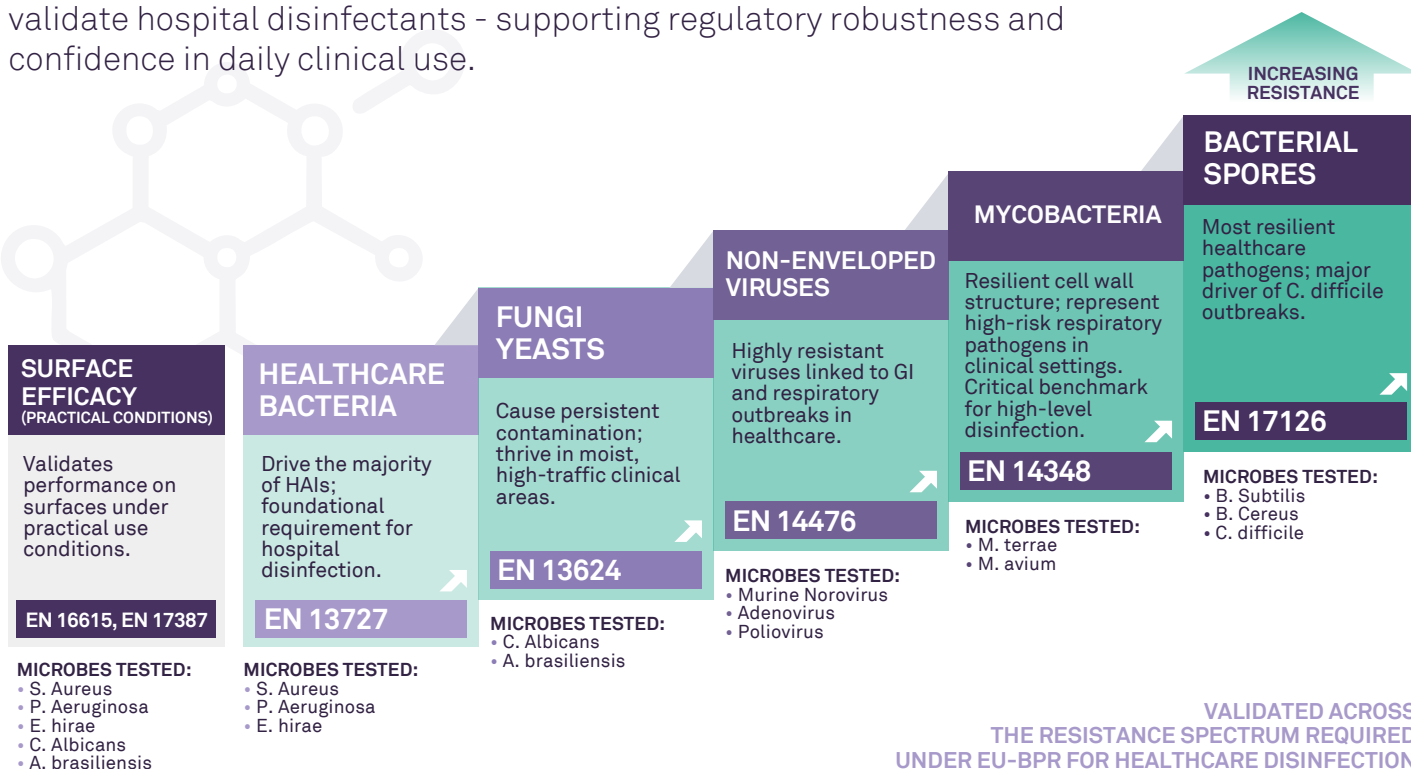

Why this matters for compliance & procurement

- Supports compliance with current EU regulatory requirements
- Reduces risk associated with outdated efficacy standards
- Provides clear evidence to support procurement evaluations
- Meets the needs of infection control, hygiene, and purchasing stakeholders



Broad-Spectrum Efficacy for Real-World Clinical Conditions

- ✓

Proven in routine healthcare settings
- ✓

Tested under high organic load
- ✓

Aligned with routine hospital workflows
- ✓

Suitable for high-risk and critical care areas

The following application guidance reflects the concentrations and contact times specified in the official Instructions for Use (IFU), aligned with EU-BPR efficacy requirements.

EN Standards	Application ²	Concentration (ppm)	Contact Time	Tablets Needed 2.5g	Tablets Needed 5g
EN 13727, EN 13624, EN 14476, EN 14348, EN 17126, EN 16615	Blood Spillage	10,000	15 mins	7 tablets in 1L water	9 tablets in 2.5L water
EN 13727, EN 13624, EN 14476, EN 14348, EN 17126	Pipette Jars	1,000	1 hr	4 tablets in 5L water	4 tablets in 10L water
EN 13727, EN 13624, EN 14476, EN 14348, EN 17126, EN 16615	General Lab/ Environmental use	1,000	1 hr	4 tablets in 5L water	4 tablets in 10L water
EN 13727, EN 13624, EN 14476, EN 14348, EN 17126	Soiled or Infected Linens	2,500	1 hr	9 tablets in 5L water	9 tablets in 10L water
EN 13727, EN 13624, EN 17387	Work Surfaces, Cupboards, Floors in Operating Theatres, ICUs, Critical Care Areas, Hallways, Bathrooms/ Washrooms	2,500	10 mins	9 tablets in 5L water	9 tablets in 10L water
EN 13727, EN 13624, EN 14476, EN 14348, EN 17126	Dish clothes, Mops, etc.	2,500	1 hr	9 tablets in 5L water	9 tablets in 10L water

PRODUCT NAME	NUMBER	PACKAGING
PRESEPT™ Tablets - 2.5g	SPT25	100 Tabs/Tub, 6 Tubs/Case
PRESEPT™ Tablets - 5g	SPT50	50 Tabs/Tub, 6 Tubs/Case

