

STEAM

INTENSIFYING OUR FOCUS ON **STERILITY ASSURANCE**

In diverse healthcare settings, the significance of steam sterilization for safety cannot be overstated. Whether in a sprawling university hospital or an ambulatory surgery center, both face mounting pressure for quicker sterile instrument turnaround, making prompt sterilization results crucial.

BIOTRACE™

AUTO READ 20 STEAM BIOLOGICAL INDICATOR
AUTO READ 20 STEAM BI PROCESS CHALLENGE DEVICE (PCD)

The **BIOTRACE™** Auto Read 20 Steam **Biological Indicator's** read-out capability, as quick as 20 minutes, minimizes downtime and allows for efficient monitoring of sterilization processes.

PRODUCT # 73100

The Auto Read 20 Steam **BI Process Challenge Device** poses a challenge to the sterilization process by demonstrating resistance equivalent to the ANSI/AAMI ST79 16 towel pack. Compliant with ISO 11140-1 (Type 2).

PRODUCT # 73130



VERISURE™

STEAM TYPE 5 MIGRATING INTEGRATORS

The **VERISURE™** Steam Type 5 Migrating Integrator combines a steam and temperature-sensitive chemical pellet with a paper wick in a laminate to monitor time, temperature, and steam, providing fast and easy-to-read results.


PRODUCT # 73150

Steam Type 5 Migrating **Integrator Process Challenge Device** poses a challenge to the sterilization process by demonstrating resistance equivalent to the ANSI/AAMI ST79 16 towel pack. Compliant with ISO 11140-1 (Type 2).

PRODUCT # 73160

Steam Type 5 Migrating Integrator with Extender for easy accessibility within surgical tray.

PRODUCT # 73155

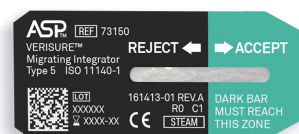

VERISURE™ BOWIE-DICK TEST PACK

VERISURE™ Bowie-Dick Test (ISO 11140-1 Type 2) Packs are designed to evaluate the efficiency of air elimination within vacuum-assisted steam sterilizers.

PRODUCT US # 73180 PRODUCT EU # 73181 *not FDA cleared



SEALSURE™ STEAM INDICATOR TAPE 18MM

SEALSURE™ Steam Indicator Tape 18mm is a steam sterilization exposure monitor which is designed to securely seal sterilization packs and provide visual assurance to SPD and OR staff to distinguish between exposed and non-exposed packs.

PRODUCT # 73170



PRODUCT #	PRODUCT DESCRIPTION	PACKAGING*
73100	BIOTRACE™ Auto Read 20 Steam Biological Indicator	50 BIs per box / 2 boxes per case
73130	BIOTRACE™ Auto Read 20 Steam BI Process Challenge Device	25 PCDs + 25 BIs (controls) per case
73150	VERISURE™ Steam Type 5 Migrating Integrator	250 CIs per bag / 2 bags per case
73160	VERISURE™ Steam Type 5 Migrating Integrator Process Challenge Device	25 PCDs per case
73155	VERISURE™ Steam Type 5 Migrating Integrator w/ Extender	100 CIs per bag / 2 bags per case
73180	VERISURE™ Bowie-Dick Test Pack US Model	20 test packs per case
73181	VERISURE™ Bowie-Dick Test Pack EU Model	20 test packs per case
73170	SEALSURE™ Steam Indicator Tape 18mm	50m per roll / 4 rolls per case

*the minimum order quantity is one case.

All CE marked products are compliant to necessary standards in order to ensure product safety and performance. Biological indicators comply with ISO 11138; chemical indicators comply with ISO 11140. FDA cleared unless otherwise stated.

asp.com