



PRESEPT™

EFFERVESCENT DISINFECTANT TABLETS



Broad-spectrum disinfection you can trust for healthcare, laboratory, and institutional hygiene applications.

ASP

PRESEPT™

EFFERVESCENT DISINFECTANT TABLETS



The Science Behind

PRESEPT™ Effervescent Disinfectant Tablets are high-level disinfection tablets based on sodium dichloroisocyanurate (NaDCC), a proven chlorine-releasing compound.

When dissolved in water, the **PRESEPT™** Solution generates hypochlorous acid (HOCl), the active antimicrobial agent responsible for rapid microorganism inactivation.

Through a unique equilibrium (reservoir) effect, **PRESEPT™** Disinfectant Tablets maintain a balance between free available chlorine and a bound chlorine reserve, allowing active chlorine to be continuously replenished as it is consumed by microorganisms or organic matter. This controlled release mechanism supports sustained antimicrobial activity and consistent performance throughout the intended contact time.

WHY NaDCC MATTERS

NaDCC is the active chemistry behind **PRESEPT™** Disinfectant Tablets, delivering controlled chlorine release and performance advantages over traditional hypochlorite-based disinfectants.

The prepared **PRESEPT™** Solution provides fast, broad-spectrum antimicrobial action at a near-neutral pH, helping reduce corrosivity while improving safety for both users and equipment.

Stabilised by monosodium cyanurate, this chemistry provides extended residual efficacy and consistent performance throughout the disinfection process.





PRESEPT™

EFFERVESCENT DISINFECTANT TABLETS



Features & Benefits

- ✓

Broad-Spectrum Efficacy
Proven activity against bacteria, fungi, viruses, mycobacteria, and spores
- ✓

Scientifically Validated Efficacy
Demonstrated performance across key pathogen groups relevant to healthcare disinfection.
- ✓

Stable Chlorine Release
Consistent potency & long shelf life for reliable disinfection
- ✓

Near-Neutral pH (6–7)
Safer for personnel and less corrosive than bleach
- ✓

Easy to use
Effervescent tablets enable simple, safe preparation with accurate dosing and rapid deployment
- ✓

Wide Range of Applications
Suitable for use across healthcare, laboratory, emergency, and institutional hygiene settings
- ✓

Multi-Year Shelf Life
Reduces waste and inventory costs with a multi-year shelf life
- ✓

Sustainable Packaging
Reduces packaging waste and storage requirements

Where PRESEPT™ Disinfectant Tablets Are Used*


Healthcare & Hospitals	Blood spill management, surface disinfection, equipment disinfection, and linen processing across hospital environments, including maternity units, nurseries, surgical areas, wards, and support areas.
Emergency Response	Rapid, reliable disinfection for blood spills, contamination events, and urgent hygiene needs. Portable tablets enable fast preparation and dependable chlorine release.
Laboratory & Research Facilities	Disinfection of glassware, pipette jars, work surfaces, and laboratory equipment. Suitable for controlled chlorine preparation in laboratory and environmental use.
Facility & Environmental Hygiene	Routine disinfection of work surfaces, cupboards, floors, bathrooms, shared areas, and cleaning tools. Also suitable for dishcloths, mops, and general environmental sanitation.
Food Preparation & Utensil Disinfection	Suitable for crockery, eating utensils, baby bottles, teats, stainless steel utensils, porcelain, glassware, and washable kitchenware in healthcare and institutional environments.





PRESEPT™

EFFERVESCENT DISINFECTANT TABLETS

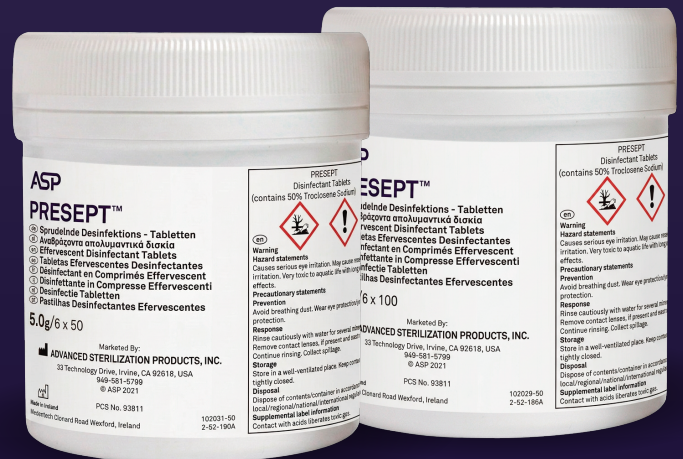

PRODUCT NAME	PRODUCT NUMBER	PACKAGING
PRESEPT™ Tablets - 2.5g	SPR25	100 Tabs/Tub, 6 Tubs/Case
PRESEPT™ Tablets - 5g	SPW50	50 Tabs/Tub, 6 Tubs/Case


FOR MORE INFORMATION VISIT

asp.com

OR CONTACT YOUR LOCAL ASP REPRESENTATIVE



ASP International GmbH, Im Majorenacker 10,
Schaffhausen Switzerland
©ASP 2026. All Rights Reserved.

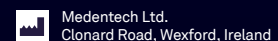

*Applications shown are aligned with regional Instructions for Use.

Prior to use, refer to the complete instructions for use supplied with the device(s) for proper use, indications, contraindications, warnings and precautions.
This document is approved for distribution in non-EU countries only.

AP-2600098-01-1